If you live with a condition that makes it harder to clear mucus, breathing can feel more difficult.

This may include:



Bronchiectasis

Cystic fibrosis

Spinal-cord injury

Neuromuscular conditions

BiWaze® Clear is designed to help loosen and move mucus in a treatment that fits more easily into your daily routine. Its dual-blower system gently opens your airways and loosens trapped mucus in a typical 10-minute session.*

By moving mucus toward the larger airways, BiWaze Clear may make it easier to cough or suction it out, so you can breathe more comfortably.

With a lightweight design and a built-in battery, you can take your treatment from the living room to the backyard — or even on the road — without missing a session.




Included with every BiWaze® Clear system:

- Control unit with an internal battery which lasts up to 10 standard therapies on a full charge*
- Convenient carrying case
- Dual Lumen Breathing Circuit™ which includes a coaxial filter, coaxial breathing tube, Aerogen® nebulizer, handset, and patient interface (mouthpiece, face mask or trach adapter)
- Power cables for the control unit and Aerogen® solo nebulizer
- Instructions for use manual and quick reference guide



Stay connected with the free Arc Connect™ app

Arc Connect™

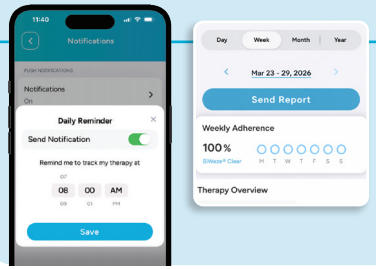

Track therapy, view progress, set reminders, and share updates.

Scan to download

RESOURCES

Website | Videos | Arc Connect™ App

* Nominal = a 10 minute therapy.
** FDA Indication to use the Aerogen nebulizer with saline solution only.

℞ Only

abm RESPIRATORY CARE

BiWaze® and ABM Respiratory Care® are registered trademarks of ABMRC, LLC.
©2026 ABMRC, LLC. ALL RIGHTS RESERVED.
PRTN-1585539264-871 Rev 3.0 SEPT-2026

For further information about this product or a service, visit our webpage: www.abmrc.com or call customer support: 1-877-ABM-RC01 (877-226-7201)

BiWaze® Clear SYSTEM

Information for Patients & Caregivers



Respiratory Therapy

3 therapies

1 system

10 minutes*

BiWaze® Clear combines three airway-clearance therapies in one treatment to help loosen and move mucus so it's easier to clear.

OPEN

Positive Expiratory Pressure (PEP)

Gentle pressure helps keep your airways open, creating a path for mucus to move.



SHAKE

High-Frequency Oscillation

Rapid pulses of air help loosen mucus and move it toward the larger airways, where it can be coughed or suctioned out.



THIN

Nebulization

Aerosol can be delivered during therapy to help thin and loosen mucus, making it easier to move and clear.



Together, these three therapies help make airway clearance simpler and more efficient **— all with one system.**

Internal or External Oscillation — Which May Be Right for You?

Your doctor may prescribe airway-oscillation therapy for reasons such as thick mucus that is hard to cough up, frequent lung infections, or limited mobility after surgery or illness.

How do the two options compare?

External oscillation (HFCWO) — vest therapy

- A vest or wrap gently vibrates your chest for about 30 minutes to loosen and move mucus.
- May be a good fit if you prefer not to use a face mask, mouthpiece, or trach adapter.

Internal oscillation (OLE) — BiWaze® Clear

- Gentle air pulses through a face mask, mouthpiece, or trach adapter help clear mucus in **about 10 minutes**, while also providing lung expansion and nebulization.
- May be a good fit if you:
  - Find a vest uncomfortable due to incisions, ports, osteoporosis, or chest sensitivity
  - Want to combine multiple airway-clearance therapies to shorten treatment time

Work with your clinician to decide which option is the right fit for your care.

Also pairs with the free **Arc Connect™** app to track therapy, view progress, set goals and reminders, and share updates if you choose.

